

BACTRACK® SingleShot™

0.08%
TESTER

BREATH ALCOHOL TEST

- Easy to Use
- Results in 2 Minutes
- Over 99% Accurate!


A change in color indicates that your Blood Alcohol Content (%BAC) is above the specified BAC level.

See top right corner for BAC test value of this kit.



Don't Drink & Drive!

Stay Safe and Consume Alcohol responsibly

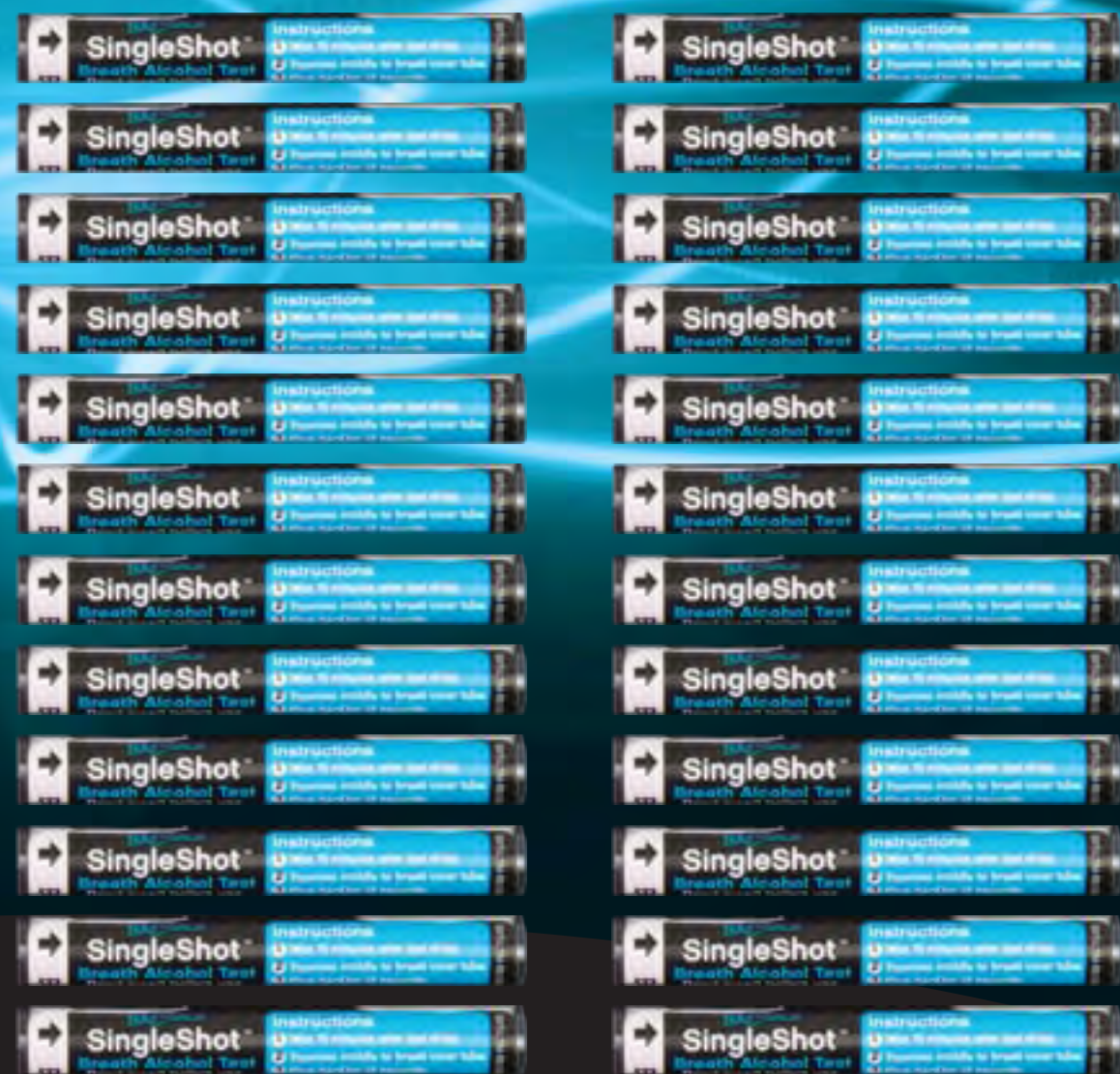

24 PACK

Great for personal use, schools,
businesses, clinics, and more!

BACTRACK® SingleShot™

Keep yourself, friends, and family safe with the BACtrack® SingleShot™. Simply blow through the tube to test whether a person's alcohol level is above the Blood Alcohol Content (%BAC) indicated on this test kit.

EASY TEST PROCEDURE*

- | | |
|---------------|--|
| STEP 1 | Wait 15 Minutes After Last Alcoholic Drink. |
| STEP 2 | Squeeze Only Once, Do Not Crush or Bend Tube. |
| STEP 3 | Hold Vertically So Arrows Are Pointing Down and Blow for 12 Seconds. |
| STEP 4 | Read The Results. |


The BACtrack SingleShot is available at 0.02, 0.05, and 0.08 %BAC thresholds. The 0.02 test is perfect for zero-tolerance environments like schools, businesses, and for family use. Higher 0.05 and 0.08 thresholds can help indicate increased levels of impairment and help prevent alcohol-related accidents.

CONTENTS

24 SingleShot™ Breath Alcohol Tubes and 1 Instruction Insert

The SingleShot is a test for alcohol in human breath. It is a disposable screening device for one-time use.

**Stay Safe and Make Smart
Decisions with the SingleShot!**

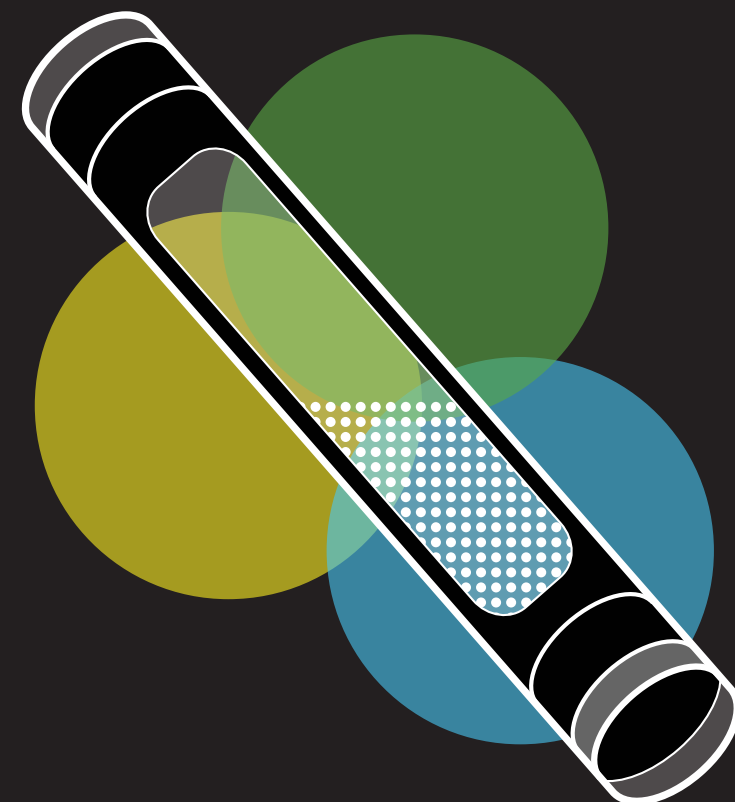

Positive Result
Green/Blue

Negative Result
Yellow

For in vitro diagnostic use.
For external use only.
Store at 4°F to 122°F
(-15°C to 50°C)
*Read Full Instructions
Prior To Use.

Distributed by:
KHN Solutions LLC
300 Broadway, Suite 26
San Francisco, CA 94133
877.334.6876
© 2009 KHN Solutions LLC



BT-SS2PK08

Lot and Expiry
Info Area